



1000-DR

Digital Conversion Package for Your PC-1000




1000-DR Digital Conversion Package

Turn your PC-1000 into a state-of-the-art direct digital panoramic X-ray system. For dental professionals already using the PC-1000, the 1000-DR Digital Conversion Package offers the quality, speed and ease that make it the right digital solution for your practice.

The 1000-DR Digital Conversion Package includes:

- Cadmium Telluride (CdTe) digital sensor.
- Dell computer system and diagnostic software.
- Same-day training (in-office or web and phone).
- Integration with your computer network, and image and practice management software.
- Single source of support, service and accountability.
- Comprehensive 2-year warranty on entire imaging system, including your existing machine.

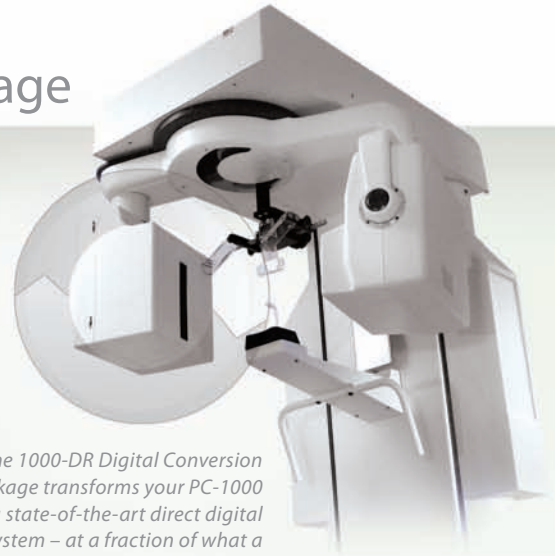

The 1000-DR Digital Conversion Package transforms your PC-1000 into a state-of-the-art direct digital system – at a fraction of what a new digital unit might cost.

Quality

- Improved diagnostics.
- CdTe technology captures images with improved clarity vs. other digital technologies and reduces the spinal column shadow typical in film X-rays.
- Diagnostic software provides multiple magnification options and further image enhancement.
- Image data allows selective enhancement of different image layers, aiding diagnosis.
- Exceptional service and technical support.
- A nationwide network of trained field service technicians to ensure your needs are promptly met.

- Free telephone and internet support from knowledgeable factory service technicians.
- Legendary durability backed by the industry's longest and most exclusive warranty.

Speed

- Conversion process, including training, takes just a few hours.
- Images are immediately available for diagnostics/patient consultation through the provided Dell® computer system with flat panel monitor.
- Patient experience and capacity are improved by a faster overall process.

- Insurance payments are often received quicker when claims are submitted electronically.

Ease

- Familiar functions and patient procedures mean your staff already has many of the skills needed, requiring a minimal learning curve.
- Integration with computer network and practice management software simplifies front office procedures and enhances patient consultation.
- Automatic image optimization provides superb images instantly with little or no need for post-capture adjustments.

Panoramic Platinum™ Care Package

With the purchase of any Digital Imaging System from Panoramic, our customers have access to one of the most comprehensive support and service packages in the industry. We call it the Platinum Care Package.



Live Phone Support



Local Technicians



Product Support



Live Web Support



Software Upgrades



Live Chat




Left: Digital Panoramic taken on a 27-year-old male.




Right: Digital Panoramic taken on a nine-year-old male.

The Cadmium Telluride (CdTe) sensor technology is more efficient because it minimizes light distortion by directly converting X-ray energy into electronic digital output. As a result, image clarity and contrast are improved up to 300%. In addition, the CdTe sensor produces over 300 image frames per second to produce a multi-layered panoramic image of unparalleled quality and clarity. This depth of image data allows selective enhancement of the different layers. As a result, poor patient positioning or unclear anatomical landmarks can easily be resolved.

Conventional CCD/CMOS

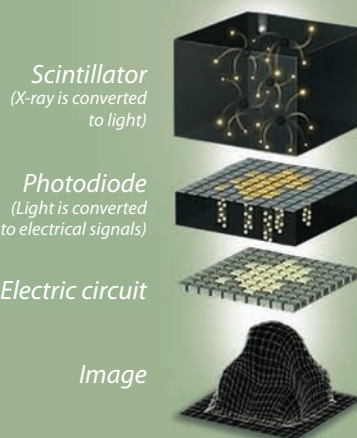

Panoramic CdTe

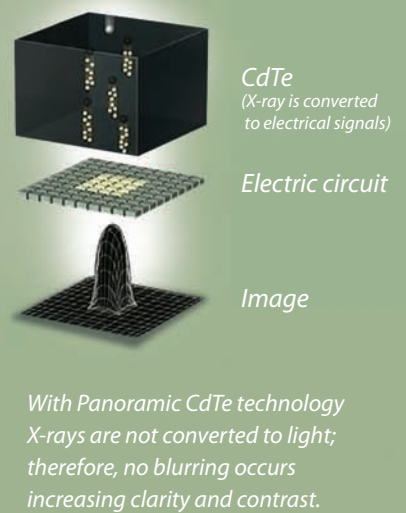

"Anytime you get a combination of great value, great service and a fantastic product, you need to celebrate! Our office truly appreciates our digital pan from Panoramic Corporation. We are getting better images than ever, with absolutely no cost to take a panorex at any time! No film to worry with, no accidental exposure of the box of film, just consistently great images. And it works great with our imaging system as well. There just aren't that many no-brainers out there any longer, but buying this unit is truly a no-brainer."

*Dr. Richard Coker - Tyler, Texas
Digital Panormaic Customer*



Panoramic Corporation
4321 Goshen Road
Fort Wayne, IN 46818
800-654-2027
www.pancorp.com